INSTRUCTIONS FOR USE

Intended Use

TransFix®/EDTA Cerebrospinal Fluid (CSF) 1-4ml Sample Storage Tubes are intended for the stabilisation and storage of human CSF specimens for characterisation of infiltrated white blood cells by flow cytometry. Recovery of infiltrated leucocyte subsets can be accomplished over a 3-day period following stabilisation.

TransFix/EDTA CSF 1-4ml Sample Storage Tubes are for **In Vitro Diagnostic Use**.

Summary and Principles

Within CSF specimens, the qualitative and quantitative characterisation of infiltrated leucocytes via flow cytometry is an important tool for identification of various immunodeficiencies and haematologic diseases [1-4]. However, the number of leucocytes present in CSF is often very low and degrade quickly. Therefore, assessment of these cells is reliant upon urgent analysis which is not always possible.

Subsets of leucocytes can be distinguished by cell surface antigens using fluorescent antibodies and flow cytometry. The TransFix stabiliser acts by preserving cell surface antigens of lymphocyte subsets until processing and analysis can be performed.

Addition of CSF to TransFix/EDTA CSF 1-4ml Sample Storage Tubes at the time of lumbar puncture has been shown to significantly extend the integrity of these leucocytes compared to unstabilised samples or stabilisation by other methods [1, see Representative Data Figure 1]. TransFix/EDTA preserves the cell surface antigens until processing and analysis can be performed [1, 2].

TransFix/EDTA CSF 1-4ml Sample Storage Tubes consist of 5ml polypropylene vials containing 0.2ml TransFix/EDTA suitable for stabilising 1-4ml of CSF.

Precautions and Warnings

1. TransFix/EDTA CSF 1-4ml Sample Storage Tube is intended for use as specified in this document. It is an 'In Vitro Diagnostic' product for professional use only.
2. TransFix/EDTA CSF 1-4ml Sample Storage Tubes are only for use with flow cytometry analysis.
3. TransFix/EDTA CSF 1-4ml Sample Storage Tubes must be stored at 2-8°C. Do not freeze.
4. Use the correct volume of CSF to ensure accurate results.
5. Do not use after the expiration date on the tubes and packaging.
6. Do not dilute or add other components to the tube before use.
7. Do not use cell viability stains on samples treated with TransFix as they are fixed instantaneously.
8. TransFix/EDTA CSF 1-4ml Sample Storage Tubes are single-use only. Do not re-use.
9. TransFix treated samples and all materials coming into contact with it should be handled as if capable of transmitting infection.
10. Avoid contact of TransFix treated samples with the skin and mucous membranes. The cell preservative is an irritant and any contact should be washed off with soap and water immediately.
11. Safety Data Sheet can be obtained at www.cytomark.co.uk or by calling +44(0)1280 827460.

Reagents

TransFix is a clear green liquid containing formaldehyde and other chemicals.

Indications of Product Deterioration

1. Cloudiness or precipitate visible in the TransFix.
2. Colour change of TransFix from a clear green liquid.
3. TransFix change from liquid to solid.

If indications of product deterioration occur, do not use and contact Cytomark immediately on: +44(0)1280 827460 or support@cytomark.co.uk.

Pre-Use Storage Conditions and Stability

TransFix/EDTA CSF 1-4ml Sample Storage Tubes are shipped in ambient conditions, between temperatures of -5°C to 26°C. Additional insulation may be required for shipment during extreme temperature conditions. TransFix/EDTA CSF 1-4ml Sample Storage Tubes must be stored at 2-8°C on arrival.

Unused TransFix/EDTA CSF 1-4ml Sample Storage Tubes are stable at 2 - 8°C for up to 12 months or until the expiration date on the label.

TransFix/EDTA CSF 1-4ml Sample Storage Tubes are supplied in sealed foil pouches.

Instructions for Use

1. Collect cerebrospinal fluid (CSF) by lumbar puncture according to the relevant clinical standard.
2. Carefully transfer 1 – 4ml CSF into the 5ml tube using a manual pipette as soon as possible and no later than 1 hour after lumbar puncture.
3. Gently mix the TransFix treated sample and close the cap.
4. Store / transport the TransFix treated sample for up to 72 hours at 2 - 8°C.
5. Stain and analyse your stabilised CSF sample as per your usual CSF flow cytometry protocol (N.B. see Notes)

Notes

1. It is recommended that all antibody conjugates are validated in association with TransFix/EDTA CSF 1-4ml Sample Storage Tubes prior to use.
2. Use caution when implementing automatic gating strategies as light scatter positions of cells stabilised by TransFix may differ from those of untreated cells.
3. Gating for lymphocytes should be performed on CD45+/SSc.
4. If calculating absolute cell counts, the dilution with TransFix/EDTA must be accounted for. I.e., for a sample size of 1ml adjust the absolute cell count by multiplying the output by 1.2. For a sample size of 2ml adjust the absolute cell count by multiplying the output by 1.1, etc.

Disposal

TransFix (EU Waste Code 18-01-06) contains formaldehyde which should be disposed of in accordance with local regulations. Avoid disposing into drainage systems and the environment.

Once the TransFix/EDTA CSF 1-4ml Sample Storage Tube contains a biological sample it must be considered an 'Absolute Hazard' (EU Waste Code 18-01-03) and disposal is in accordance with local regulations regarding clinical waste.

References

1. Use of TransFix Cerebrospinal Fluid Storage Tubes Prevents Cellular Loss and Enhances Flow Cytometric Detection of Malignant Haematological Cells after 18 Hours of Storage. De Jongste *et. al.*, Cytometry Part B 2014; 86B: 272– 279.
2. Guidelines on the use of Multicolour Flow Cytometry in the Diagnosis of Haematological Neoplasms. Johansson *et. al.*, British Journal of Haematology, 2014, 165, 455-488.
3. Flow Cytometric Characterization of Cerebrospinal Fluid Cells. De Graaf MT *et. al.*, Cytometry, 2011, 80B: 271-281.
4. ESCCA/ISCCA Protocol for the Analysis of Cerebrospinal Fluid by Multiparametric Flow-cytometry in Hematological Malignancies. Del Principe *et. al.*, Cytometry, 2020; 1– 13.

Ordering Information


Please call Cytomark on +44(0)1280 827460 or email support@cytomark.co.uk for assistance. Additional information can be found online at www.cytomark.co.uk. A certificate of conformity can be provided with every batch of TransFix/EDTA CSF 1-4ml Sample Storage Tubes.

Product Description	Catalogue Numbers
TransFix / EDTA CSF 1-4ml Sample Storage Tube (2 tubes)	TF-CSF-L-2
TransFix / EDTA CSF 1-4ml Sample Storage Tube (10 tubes)	TF-CSF-L-10
TransFix / EDTA CSF 1-4ml Sample Storage Tube (25 tubes)	TF-CSF-L-25
TransFix / EDTA CSF 1-4ml Sample Storage Tube (50 tubes)	TF-CSF-L-50

Cytomark, a Division of Caltag Medsystems Ltd.

Caltag Medsystems Ltd, Whiteleaf Business Centre, 11 Little Balmer, Buckingham, MK18 1TF, UK

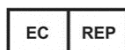Advena Ltd, Tower Business Centre, 2nd Flr, Tower Street, Swatar, BKR 4013, Malta

REPRESENTATIVE DATA (P-CSF-005 REPORT Length of Stabilisation of CSF Samples treated with TransFix/EDTA (MC009), Caltag Medsystems Ltd 2021–unpublished)

Figure 1 a-e. shows recovery of monocytes, T cell subsets, NK cells and B cells from Mock CSF (Thermo Fisher Scientific MAS-CSF (CSF-MP) spiked with PBMCs) stabilised with TransFix/EDTA, basal RPMI and complete RPMI, and compared to unstabilised Mock CSF (control) at various time points. Data from three replicates for each sample were averaged at each time point. Each replicate was 100µl of sample at approximately 1x10⁶ PBMCs/ml, so a total of approximately 100,000 PBMCs per replicate in each tube.

Figure 1 a. CD45+CD3+CD4+ T Cells

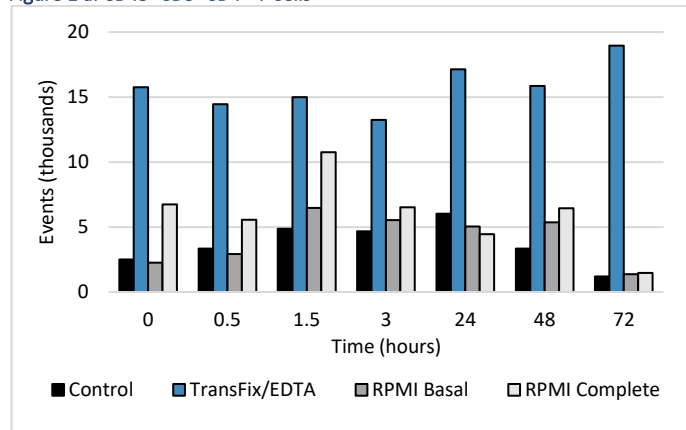
Figure 1 b. CD45+CD3+CD8+ T Cells

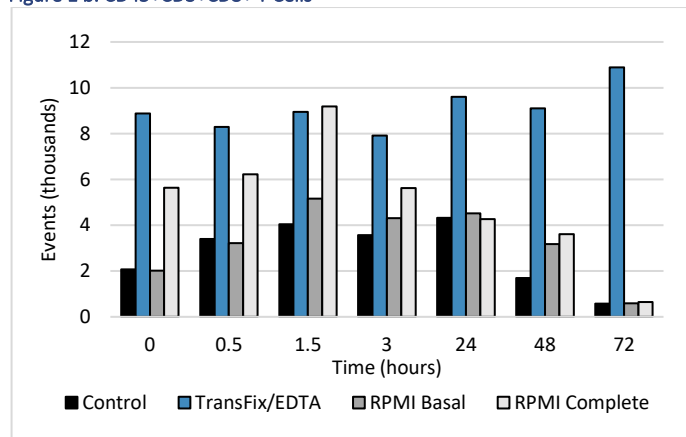
Figure 1 c. CD45+CD3-CD19+CD20+ B Cells

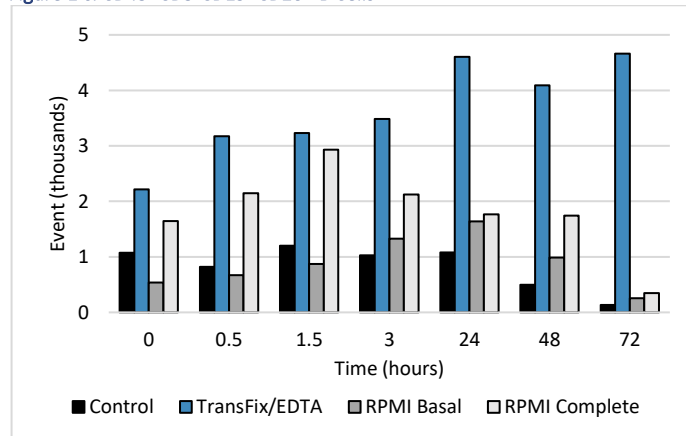
Figure 1 d. CD45+CD3-CD56+ NK Cells


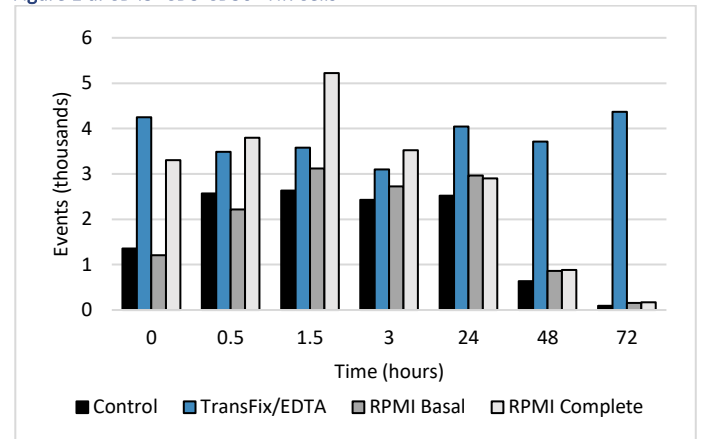
Figure 1 e. CD45+CD14+ Monocytes

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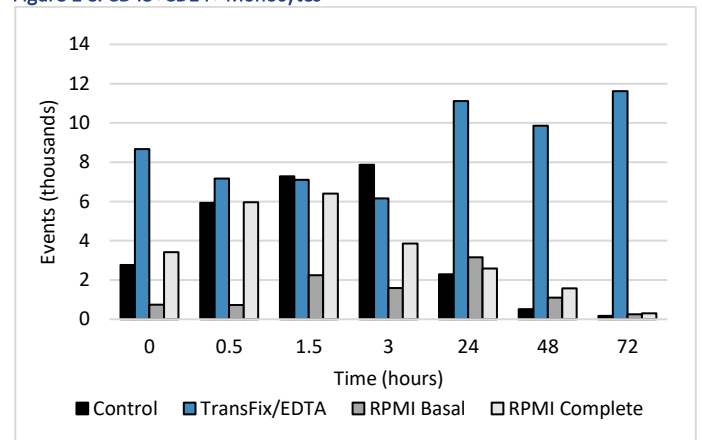Odkazy

- Use of TransFix Cerebrospinal Fluid Storage Tubes Prevents Cellular Loss and Enhances Flow Cytometric Detection of Malignant Haematological Cells after 18 Hours of Storage. De Jongste *et. al.*, Cytometry Part B 2014; 86B: 272– 279.
- Guidelines on the use of Multicolour Flow Cytometry in the Diagnosis of Haematological Neoplasms. Johansson *et. al.*, British Journal of Haematology, 2014, 165, 455-488.
- Flow Cytometric Characterization of Cerebrospinal Fluid Cells. De Graaf MT *et. al.*, Cytometry, 2011, 80B: 271-281.
- ESCCA/ISCCA Protocol for the Analysis of Cerebrospinal Fluid by Multiparametric Flow-cytometry in Hematological Malignancies. Del Principe *et. al.*, Cytometry, 2020; 1– 13.

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Riferimenti

1. Use of TransFix Cerebrospinal Fluid Storage Tubes Prevents Cellular Loss and Enhances Flow Cytometric Detection of Malignant Haematological Cells after 18 Hours of Storage. De Jongste *et. al.*, *Cytometry Part B* 2014; 86B: 272- 279.
2. Guidelines on the use of Multicolour Flow Cytometry in the Diagnosis of Haematological Neoplasms. Johansson *et. al.*, *British Journal of Haematology*, 2014, 165, 455-488.
3. Flow Cytometric Characterization of Cerebrospinal Fluid Cells. De Graaf MT *et. al.*, *Cytometry*, 2011, 80B: 271-281.
4. ESCCA/ISCCA Protocol for the Analysis of Cerebrospinal Fluid by Multiparametric Flow-cytometry in Hematological Malignancies. Del Principe *et. al.*, *Cytometry*, 2020; 1- 13.

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Caltag Medsystems Ltd, Whiteleaf Business Centre, 11 Little Balmer, Buckingham, MK18 1TF, Regno Unito

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Cytomark, en avdelning inom Caltag Medsystems Ltd.
Caltag Medsystems Ltd, Whiteleaf Business Centre, 11 Little Balmer, Buckingham, MK18 1TF, Storbritannien